



Drop-cell connector for
thin-film electrodes


Drop-cell connector



The new **drop-cell connector** (*Ref. ED-DROP-CELL*) provides a true **user-friendly** interface with the potentiostat for the thin-film (micro)electrodes supplied by MicruX Technologies.

» Drop-cell connector features

The drop-cell connector enables the use of **microvolume** ($1 - 10 \mu\text{L}$ sample drops) with all **standard** ($10 \times 6 \text{ mm}$) **thin-film (micro)electrodes** supplied by MicruX Technologies.

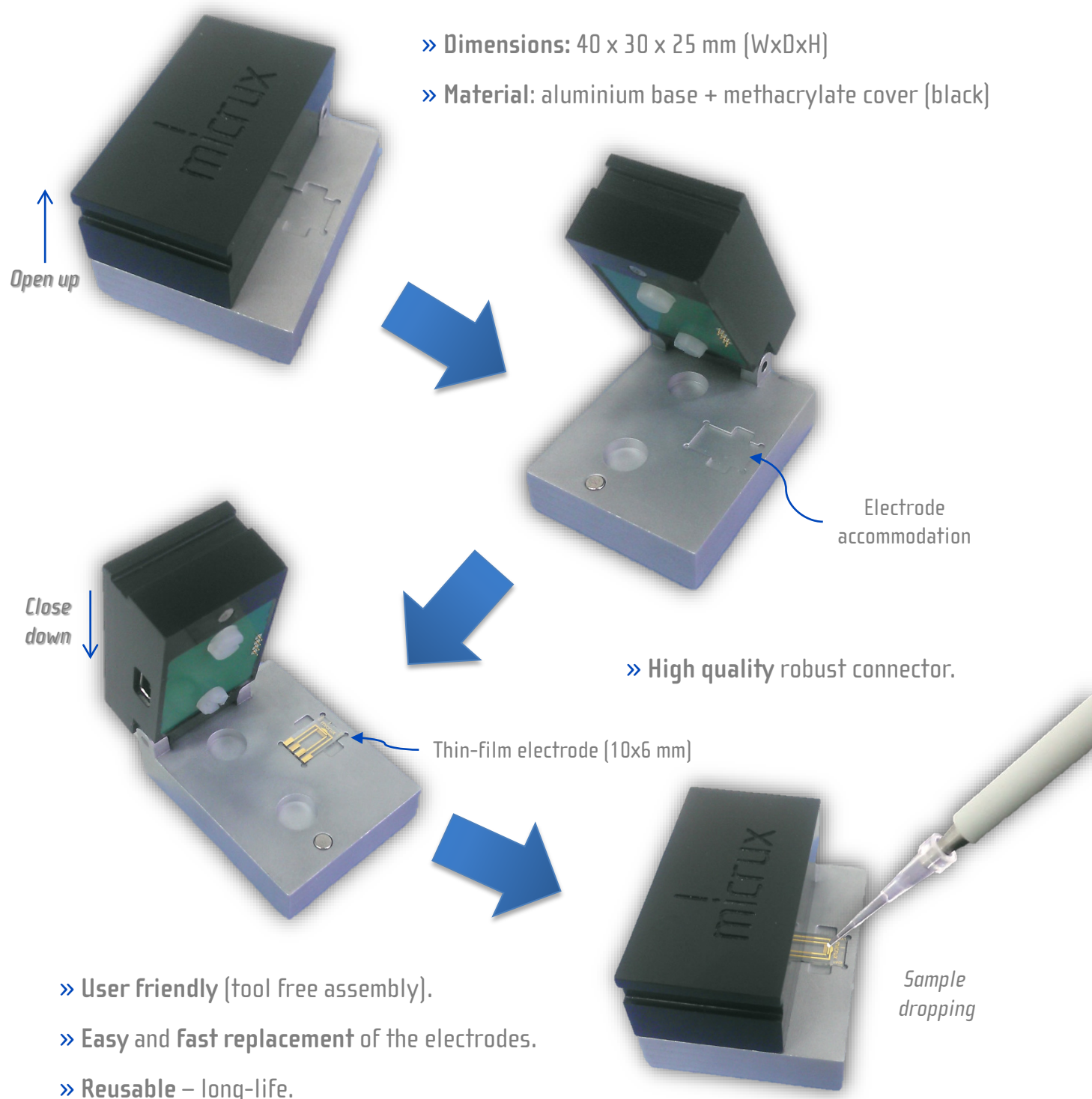

Drop-cell connector




» Drop-cell potentiostat-interfacing

» Connector back view



Integrated easy-to-use connector (miniUSB-type connector) with pogo-pins.

» Universal cable

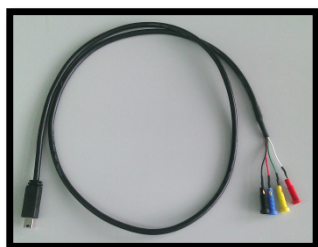


The drop-cell connector is supplied with an universal cable compatible with any commercial potentiostat.



miniUSB to EC Drop/Flow-cell

Female/male banana plug to potentiostat

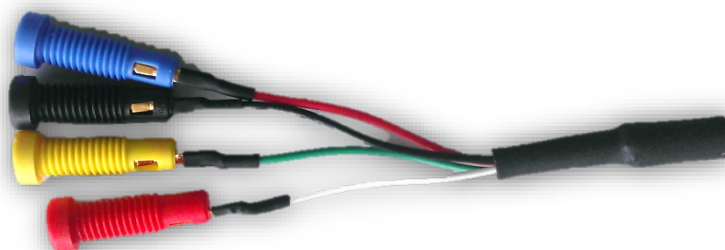

Plug to potentiostat may be available in other format under previous request.



Mora-Garay Industrial Park · Juan de la Cierva, 2C, Bldg. # 6
33211 · Gijón (Asturias) · SPAIN

Phone/FAX: +34 984151019

E-mail: info@micruxfluidic.com

Web: www.micruxfluidic.com

